

Reveal™ distal attachment cap



Use of an endoscopic cap improves visualization by “maintaining the optimal field of view and keeping the correct depth of field; pushing aside folds, angulations, tissues, or sphincters; and stabilizing the position of the tip of the endoscope.”*


**The endoscopic cap that can, Gastrointestinal Endoscopy, Volume 76, Issue 1, 172 July 2012*

Maintain a clear view during endoscopic procedures...

The Reveal™ distal attachment cap is designed to aid the endoscope during a variety of therapeutic and endoscopic procedures.

The Reveal™ distal attachment cap...

- improves visualization and helps avoid 'red out' by maintaining a constant distance between the endoscope and mucosa
- aids in flattening folds and maneuvering through difficult turns in the GI tract
- can be used in conjunction with accessories such as snares or needles

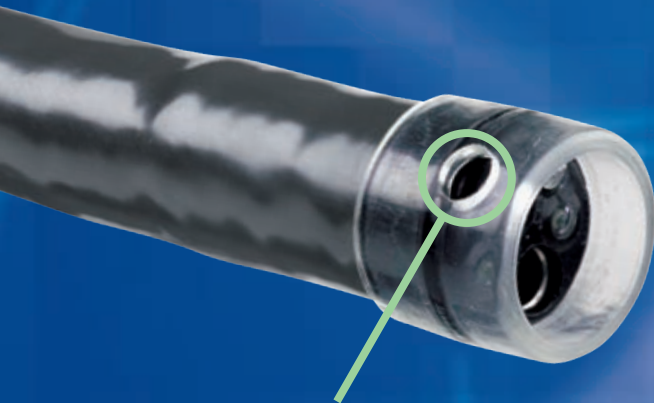

side hole allows for constant drainage of fluid

For more information on the Reveal™ distal attachment cap, contact US Endoscopy today.

product number	description	cap outer diameter (mm)	scope tip outer diameter* (mm)	length (mm)	sterile	units/box
00711771	Reveal™ distal attachment cap	11.8	9.5-9.9	4	yes	10
00711773	Reveal™ distal attachment cap	13.4	10.9-11.7	4	yes	10

*Some scope tips may have variable outer diameter tolerances.

Distributed by:



*listening...and delivering solutions**

5976 Heisley Road
Mentor, OH 44060 USA

phone +1 440 / 639.4494

fax +1 440 / 639.4495

email global@usendoscopy.com

www.usendoscopy.com

760989 Rev. A

US Endoscopy is a subsidiary of STERIS Corporation.
Listening...and delivering solutions®, the US Endoscopy listening...and delivering solutions design, and all marks denoted with ® or ™ are registered with the U.S. Patent and Trademark Office, or trademarks owned by US Endoscopy

Follow us:    